



MEDICAL DEVICES & MORE

INTENSIVE CARE PRODUCTS

LINES | PORTS



VISION

To ensure easy reach and regular availability of our products to all our customer, at all time to continue building trust among Stake holders of **SCEPTRE**.

MISSION

To serve Medical fraternity with products of international standards and create opportunities for healthy and consistent growth of our Associates, **SCEPTRE** has never lost its vision to identify the needs of Patients and serve them with most suited products.

QUALITY POLICY

By incorporating latest manufacturing techniques we at **SCEPTRE** are committed to produce best quality products with suitable prices. Fulfilment of Customer needs and exceed their expectation are our utmost priorities. This is ensured by providing environment for hassle-free work culture and uninterrupted supply of materials, components and other resources. Our goal is to further enhance markets in exports, domestic & institutional areas and improve our quality management system.



your
kind of
care

"Catheters and other key components are developed based on research done in USA, TAIWAN & PRC"

INTRODUCING NEW PRODUCTS

Spinal Needle

Endotracheal Tube

Inflation Device

Tracheostomy Tube

Alcohol Swab

Gauze Swab

Disposable Biopsy Instrument

Innovative Handle Design

Handle design is smart & functional. It is light weight for easier manipulation during procedure, provided with ribs for improved control & non-roll properties when not in use.

Instrument design is new & simpler, improved for ease of biopsy procedure with exceptional control. It offers ease of operation & enhanced control for one handed cocking. Upper slide withdraws and locks the needle cannula. The bottom slide engages and cocks the stylet.

Disposable Biopsy Instrument

Product Code	Gauge Size	Colour Sample	Needle Length	Length of Notch	Penetration Depth	Dimensions Box
301.J0001A	14 G	Green	10.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0001C	14 G	Green	16.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0002A	16 G	Purple	10.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0002C	16 G	Purple	16.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0003A	18 G	Pink	10.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0003C	18 G	Pink	16.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0003D	18 G	Pink	20.0CM	18.0 mm	22.0 mm	475x170x70mm
301.J0003E	18 G	Pink	25.0CM	18.0 mm	22.0 mm	475x170x70mm

Innovative Design • Precise • Cost Effective

Core Needle is manufactured with highest quality standard for consistent core samples with left distortion of the cellular pattern.

Minimal Tissue Displacement due to high action speed.

Cutting Cannula shears tissue sharply, capturing quality samples in needle notch.

Curved Needle Tip is grit-blasted to create more diffuse surface area.

Dual Firing option for individual preference and/or patient position.

Colour Coded Buttons for easy and accurate needle gauge identification.

Angled sample notch for enhanced needle action.



Bone Marrow Biopsy Needle

Unique Cannula Design provides outstanding quality specimen & provides an exceptional ability to capture and retain bone marrow specimens. Thus elimination painful deflection technique to harvest a diagnostic quality specimen.

Twist-Lock Handle is designed ergonomically to provide greater clinician comfort.

Stylet Point has double-diamond for easy penetration of hard bone.




Bone Marrow Biopsy Needle are available with or without a marrow extraction cannula. The former is aimed at ensuring a quality specimen by collecting & holding the bone marrow in the needle cannula without using the painful needle deflection technique or altering the sample's architecture. This increases patient comfort, reduces the risk of a lost sample, and improves the quality of the specimen.

Product Code	Size
301.H0001B	8G x 10cm.
301.H0002B	10G x 10cm.
301.H0003C	11G x 15cm.
301.H0004C	13G x 15cm.

Pressure Monitoring Line

Made of clear, non-toxic PVC.
Highly kink resistant tubing.
Facilitates optimum dynamic response.
Connector: Female, Male.
Various lengths is available.
Latex Free, Sterile, non-pyrogenic.
External Dia.: 2.80mm.
Internal Dia.: 1.27mm.
Length: 10, to 250cm.

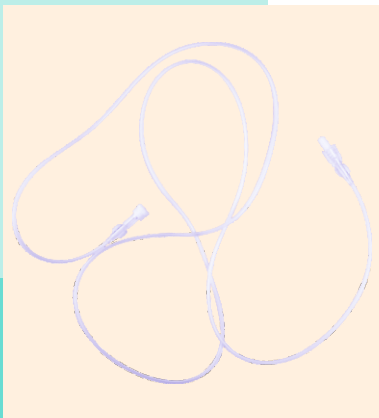

Available in Brand



Pressure Monitoring Line is made of clear and non toxic PVC. It is available as Sterile & Non-Toxic.

Tubing Material i.e., PVC is latex free & kink-resistant. It facilitates optimum dynamic response.

The Assembly is highly accurate and monitors the fluent blood pressure as required. Usually for high-risk patients the faulty readings of the instruments can be harmful. It is important in such cases that correct and timely monitoring is done for further line of treatment based on the previous readings.

Product Code	Size	Standard Packing Box/Carton	Dimentions Box/Carton
301.K0100A	10-250Cm	25/250	170x135x150mm/695x330x180mm

Intravenous Cannula

Flexible Wings are designed with 6° angle.

Stainless-steel Needle is triple facetted back cut & highly polished. It is covered with protective cap.

Catheter is kink-resistant with colour coded injection port.

Flashback Chamber is clear and transparent to ensure timely withdrawal of needle.



Product Code	Gauge & Length	Catheter O.D./Length	Standard Packing Box/Carton	Dimensions Box/Carton	Gross Wt.
301.W010HE	(14G x 1 ^{3/4"})	2.10/45	100/1000	170x135x150mm/705x425x150mm	6.50 Kg
301.W010GE	(16G x 1 ^{3/4"})	1.70/45	100/1000	170x135x150mm/705x425x150mm	6.45 Kg
301.W010FE	(17G x 1 ^{3/4"})	1.50/45	100/1000	170x135x150mm/705x425x150mm	6.40 Kg
301.W010EE	(18G x 1 ^{3/4"})	1.30/45	100/1000	170x135x150mm/705x425x150mm	6.35 Kg
301.W010DD	(20G x 1 ^{3/4"})	1.10/32	100/1000	170x135x150mm/705x425x150mm	6.30 Kg
301.W010CC	(22G x 1")	0.90/25	100/1000	170x135x150mm/705x425x150mm	6.25 Kg
301.W010BB	(24G x 3/4")	0.70/19	100/1000	170x135x150mm/705x425x150mm	6.20 Kg
301.W010AA	(26G x 5/8")	0.60/16	100/1000	170x135x150mm/705x425x150mm	6.10 Kg

Three Way Stop Cock

Tap-rotation is soft & smooth.

Uniform Cavity-bore throughout ensures maximum infusion-flow and minimum turbulence.

Leak-proof Body withstands pressure up-to ten bars when applied in normal conditions.

Luer-lock Screw ensures safe and secure connection to luers of other standard systems 6% luer-taper is compatible with any standard product.



Product Code	Color	Standard Packing Box/Carton
301.N0310A	Blue	50/1000
301.N0110A	White	50/1000
301.N0210A	Red	50/1000

Three Way Stop Cock With & Without Extension Tube



Product Code	Color	Size
301.P3120A	Blue	10 cm. - 200 cm.
301.P1120A	White	10 cm. - 200 cm.
301.P2120A	Red	10 cm. - 200 cm.

Leak-proof Three way stop cock integrated with kink resistant PVC tubing is a perfect module which facilitates drug administration or joining additional lines together. Available in various tube lengths of 10cms.-250cms.



Product Code	Size
301.P0220A	10 cm.
301.P0220C	50cm.
301.P0220E	100cm.
301.P0220G	200cm.

Kink resistant PVC tubing fitted with male luer lock at one end and female luer lock at the other is best suited to separate administration from insertion site. Available in various tube lengths of 10cms. -250cms.

Cautery Pencil

(Hand Control Pencil)

Excellent push button
disposable clearslit.
(We also have tip cleaner)
170mm - ABS pencil
appearance

Description

- Stainless Steel tip - 4cm
- Diameter of the tip - 2.3mm
- PVC medical cable - 3m
- Maximal output current - 1A
- Frequency confine - 0-1.0MHZ
- Ways of aseptic technique - Sterile EO
- Current resistance - 200 (Ohms)
- Appropriate switch sensitivity.
- Clear & appropriate switch position, with gentle "click" when activated
- Streamlined
- Anti-slippery surface
- Light weight

Perfusion Set

Roller Clamps are specially designed for better flow control & precise adjustment of infusion rate.

Spike is of super strong bevel shape and extremely sharp (Vented and Non-Vented) for better penetration.

Latex Tube is of superior quality for self-sealing.

Sharp Needle for better penetration.

Drip Chamber is designed as collapsible cylinder with hydrophobic filter that acts as bacteria barrier and is transparent to visualize the flow rate.

Options Available: Y-Injection Site / Latex-bulb / Luer-slip (with or without Air-vent)

Packing: Pouch Packed / Paper Blister Packed / DEHP Free.



Product Code	Type
301.R0110A	Vented Latex Bulb
301.R0230A	Non-Vented Y-Site

Blood Transfusion Set

Blood Transfusion Set is specially designed for transfusion of blood and blood components.

Roller Clamps are specially designed for better flow control & precise adjustment of infusion rate.

Latex Tube is of superior quality for self-sealing.

Sharp Needle for better penetration.

Drip Chamber - Single Drip / Double Drip with 200-micron filter mesh.

Options Available: Y-Injection Site / Latex-bulb / Luer-slip (with or without Air-vent)

Packing: Pouch Packed / Paper Blister Packed / DEHP Free.




Product Code	Type
301.S0110A	Vented Latex Bulb
301.S0220A	Non-Vented Y-Site

Arterial Catheter

Single Lumen Catheter

Size	Length
2.0Fr (22G)	5cm-80cm
3.0Fr (19G)	5cm-80cm
4.0Fr (18G)	5cm-80cm
5.0Fr (16G)	5cm-80cm
6.0Fr (14G)	5cm-80cm
7.0Fr (13G)	5cm-80cm
8.0Fr (12G)	5cm-80cm

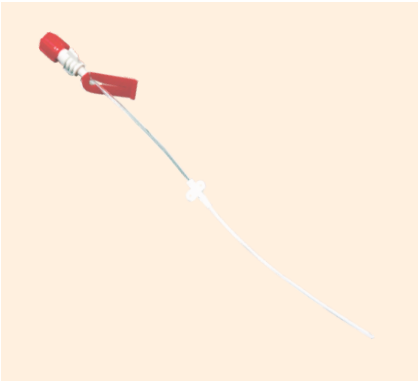

Application

- Monitor invasive arterial blood pressure
- Collect blood sample
- Infusion drug and solution

Characteristics

- The tube is good biocompatibility and it is convenient for catheter placement due to stiff catheter
- Standard seldinger technique
- Catheter and guidewire are protected by sheath
- Nitinol guidewire is available
- Straight introducer needle configuration is available

Foley Balloon Catheter


A **Foley balloon catheter** is a thin, sterile tube inserted into the bladder to drain urine. Because it can be left in place in the bladder for a period of time, it is also called an indwelling catheter. It is held in place with a balloon at the end, which is filled with sterile water to prevent the catheter from being removed from the bladder. The urine drains through the catheter tube into a bag, which is emptied when full. The procedure to insert a catheter is called catheterization.

Siliconized Latex Foley Catheter Sizes

- Foley Catheter - 6.0 Fr - 3ml
- Foley Catheter - 8.0 Fr - 3ml
- Foley Catheter - 10.0 Fr - 3-5ml
- Foley Catheter - 12.0 Fr - 5-10ml
- Foley Catheter - 14.0 Fr - 15ml
- Foley Catheter - 16.0 Fr - 30ml
- Foley Catheter - 18.0 Fr - 30ml
- Foley Catheter - 20.0 Fr - 30ml
- Foley Catheter - 22.0 Fr - 30ml

- Foley Catheter (3way)- 18.0 Fr - 30ml
- Foley Catheter (3way)- 20.0 Fr - 30ml
- Foley Catheter (3way)- 22.0 Fr - 30ml
- Foley Catheter (3way)- 24.0 Fr - 30ml

Product Code	Size
301.Z111AA	6Fr-3ml
301.Z122BB	8Fr-18ml

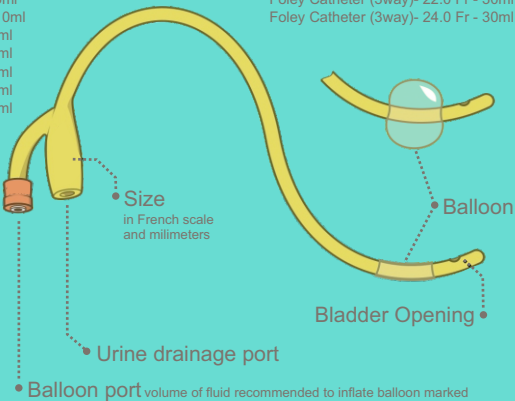

Surgical Staples

- Staples
- Staple Remover



Skin stapler is single use device.

It is designed to delevers precise rectangular shape stainless steel staple for routine wound close.

It is typically used at ent of surgical product.

It has 35 pins which is medical grade 316L SS and forms precise rectangular shape staples formation.

Implantable Subcutaneous Infusion Ports



- REDUCED WIDTH • ERGONOMIC SHAPE • LOW PROFILE
- IMPROVED AESTHETICS • MINIMUM INCISION SIZE • X-RAY VISIBLE PRINTING

Suitable for both CECT Injection and Infusion therapy.

CT Injection for diagnostic imaging (at up to 5cc/sec at 300psi).

Product Code	Size
301.4200DZ	6.6F
301.4200FZ	8.0F
301.4200HZ	9.0F

Peripherally Inserted Central Catheter

- MRI COMPATIBLE • ALCOHOL COMPATIBLE • IODINE COMPATIBLE

PICC Line are approved for both contrast / enhanced computed tomography injections and infusion therapy.

Improved Design allows for CT injections for diagnostic imaging at up to 5cc/sec at 300psi.

Available in single, dual, and triple lumen configurations.

Product Code	Size
301.B211C1	4.0Fr x 55cm (Single Lumen)
301.B211E2	5.0Fr x 60cm (Single Lumen)
301.B221E1	5.0Fr x 55cm (Double Lumen)
301.B231G2	6.0Fr x 60cm (Triple Lumen)





www.sceptremedica.com

Clinically Tested

ISO 9001 : 2015 & 13485 Certified

FDA Registered

Sceptre Medical Devices Pvt Ltd

Office : 6-Central Market, Ashok Vihar, PH-I, Delhi-110 052, DL, India

Works : 187, Sector-38, PH-I, HSIIDC, Rai, Sonapat-131 029, HR, India

Email : info@sceptremedica.com | **Tel :** +91 130 2368187

USA : 2100 North Line Street, J101 Lansdale, Pa 19446

Email : medical@prem4biz.com | **Tel :** +1(208)-861-6488

